

SAATicare

Precision Fabrics & Components for Medical & Diagnostic Applications

FOR THE MANY, BY THE FEW

—SAATI


Contents

 Click a page number to visit the section.
 To return click the house icon at the bottom of pages.

Company Information	02
Customer Service	04
SAATlcare Manufacturing	06
Fabrication Capabilities	08
Medical Applications	20
Blood Filtration	22
SAATlcare Hyphil Treatment	28
SAATlcare High Flow Fabrics	30
Diagnostic Applications	32
Testing Compliance	38
Sustainable Future	40
Interactive Product Index	42
SAATI Global Contact Information	Back

Engage
with
SAATI

SAATI S.p.A. Via Milano, 14-22070 Appiano Gentile (CO), Italy | www.saati.com

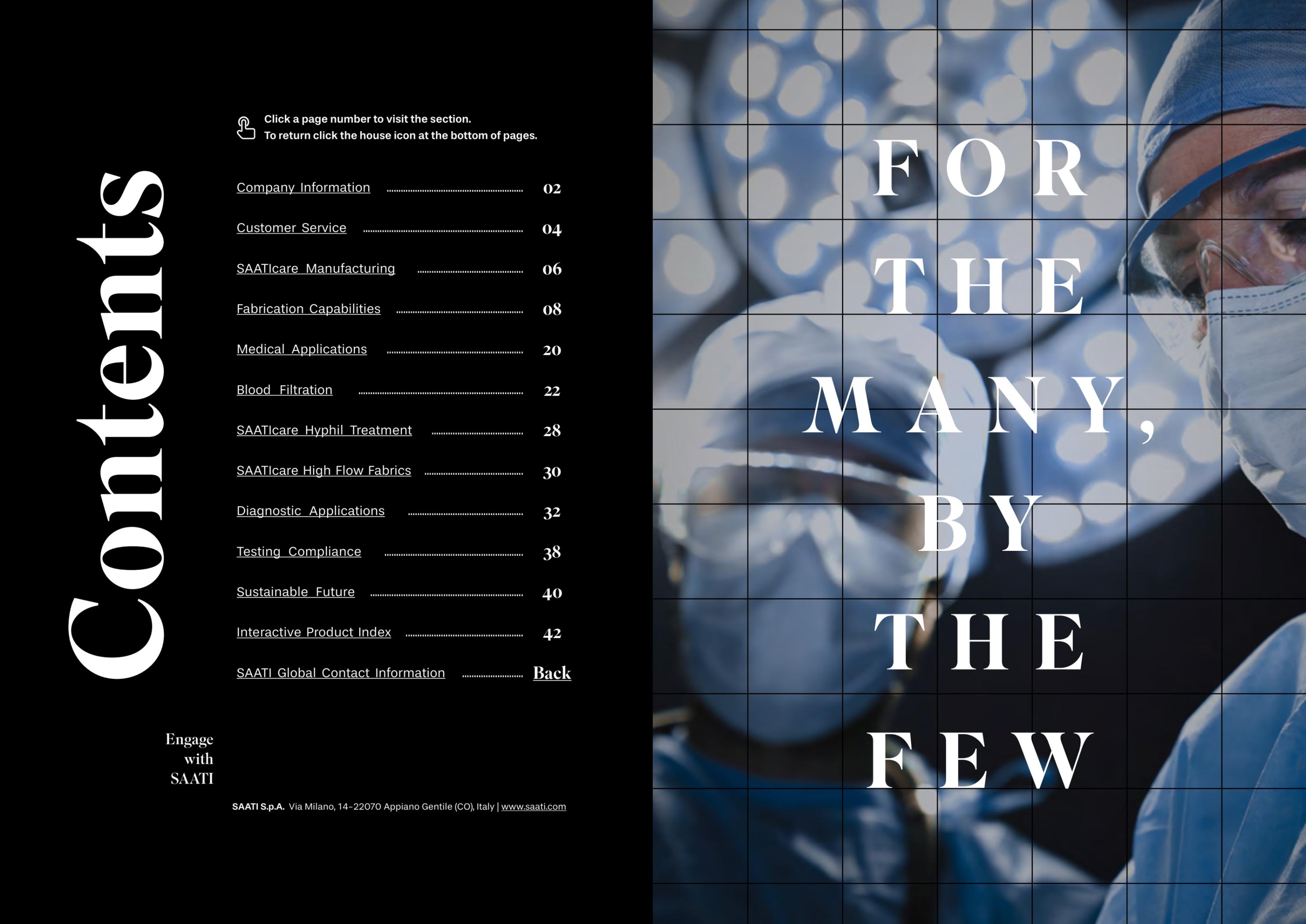

FOR
THE
MANY,
BY
THE
FEW

SAATI produces many products for all types of markets that find their way into the daily life of billions of people. We are leaders in both process and microfiltration, screen printing, and personal and vehicular multi-threat protection products.

Over Ninety Years of Innovative Action

SAATI is a multinational group with corporate headquarters that have been situated in northern Italy since 1935. Today SAATI is a leading force in the development, manufacturing and commercialization of high tech filter media & chemicals.

SAATI's passion and creativity are the foundation for an unsurpassed tradition of continuous innovation in the filtration markets. This endless pursuit is what drives SAATI's dedicated customer-centric R&D to functionalize products beyond simple filtration.

SAATI's wide range of synthetic textiles and fabricated parts in Polyester, Nylon, Polypropylene, PEEK and PPS are the ideal engineered solution for demanding microfiltration applications.

Through specialized processing and rigorous inspection, SAATI ensures consistent lot quality across tolerances, uniformity, strength, stability, and cleanliness for applications with pore sizes from 7 to 3,000 microns.

Perfecting the Art of Precision Woven Fabrics

With over 1,000 employees spread across multiple facilities worldwide, and a strong established track record in innovation and manufacturing excellence, SAATI's mission is to improve the life of every person every day - through working with both customers and partners to create a safer, healthier and cleaner world.

To guarantee product reliability, SAATI constantly runs tests and has all of the strictest and most up to date certifications that validate the consistency, performance, quality and characteristics of each item.

Specific SAATI fabrics are tested and certified in accordance with USP CLASS VI/ISO 10993 Regulations, and these fabrics are inspected and transformed into customized shapes in Class 10,000/ISO Class 7 Clean Rooms in accordance with UNI ISO 9001 regulations.



Customer Driven Innovation

Thanks to a direct presence in many countries, it is easy for customers to reach out, wherever they are located, and the response is always prompt. SAATI staff has a high level of technical expertise and dedication, and is always aiming to find the best solution for the customer's requirements.

SAATI sales representatives and engineers understand their customers' applications, working closely with staff in both the production and R&D departments to offer the customized solutions that best meets their needs.

The quality of SAATI's products is backed by the dedication and expertise of SAATI's customer service. Thanks to offices, warehouses, storage and fabrication facilities throughout the world, SAATI provides strong local support, expert responses to customer inquiries, strong engineering capability, technical support and fast delivery around the world.

SAATCicare Manufacturing

Medical and Diagnostic Devices

Medical and diagnostic devices are among the most demanding fabric applications.

SAATCicare healthcare fabrics are more than up to the task; they are the preferred choice for medical devices such as infusion and transfusion filter, arterial filter, oxygenators, micro-aggregate filters, blood bags and diagnostic applications such as test strips, spirometers, biopsy bags and molecular sieves.


The precise aperture size, uniform high flow rates, and lot-to-lot consistency make SAATCicare fabrics the ideal solution for demanding healthcare applications.

Manufactured To The Most Exacting Standards

SAATCicare fabrics are woven with monofilament polyester or polyamide fibers with smooth and uniform surfaces that are particularly suitable for medical applications. The monofilament fibers are non-shedding to reduce the risk of particulate contamination. The fibers are woven to exceptionally tight tolerances, creating uniform pore sizes, excellent strength, and good dimensional stability.

Customized to Meet Specific Customer Needs

Typical surface modifications requested include plasma treatment, adding hydrophilic or hydrophobic characteristics, and dyeing in virtually any color. SAATI Engineers can also evaluate developing a fully customized solution using your specific chemistry applied to the mesh surface to functionalize the surface beyond simple filtration action.

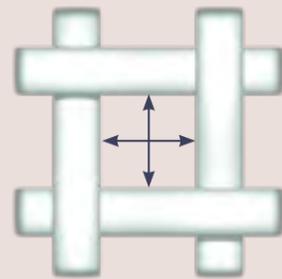

Precise Mesh Opening

The Mesh Opening is the square space between two warp and two weft yarns



High flow rates with low pressure loss



**SAATI Engineers
can also evaluate
the development of
fully customized
solutions using
your specific needs
to functionalize
the surface
beyond typical
filtration action.**


SAATCcare Fabrication Capabilities

High-Quality, Cost-Effective
Custom Fabrication from SAATI

SAATCcare fabrics can be delivered in a wide variety of finished or partly finished products. Specially designed processing equipment creates quality custom-fabricated parts in the most cost-effective way. SAATCcare fabrics can be cut into ribbons, formed into tubes, pleated, or made into any required shape. Critical fabricated parts are undergo manufacture in SAATI's Class 10,000 Clean Rooms (ISO 7).

When SAATI engineers design a part for custom fabrication, their primary consideration is to produce a high-quality precision product that will remain consistent from order to order, year after year.




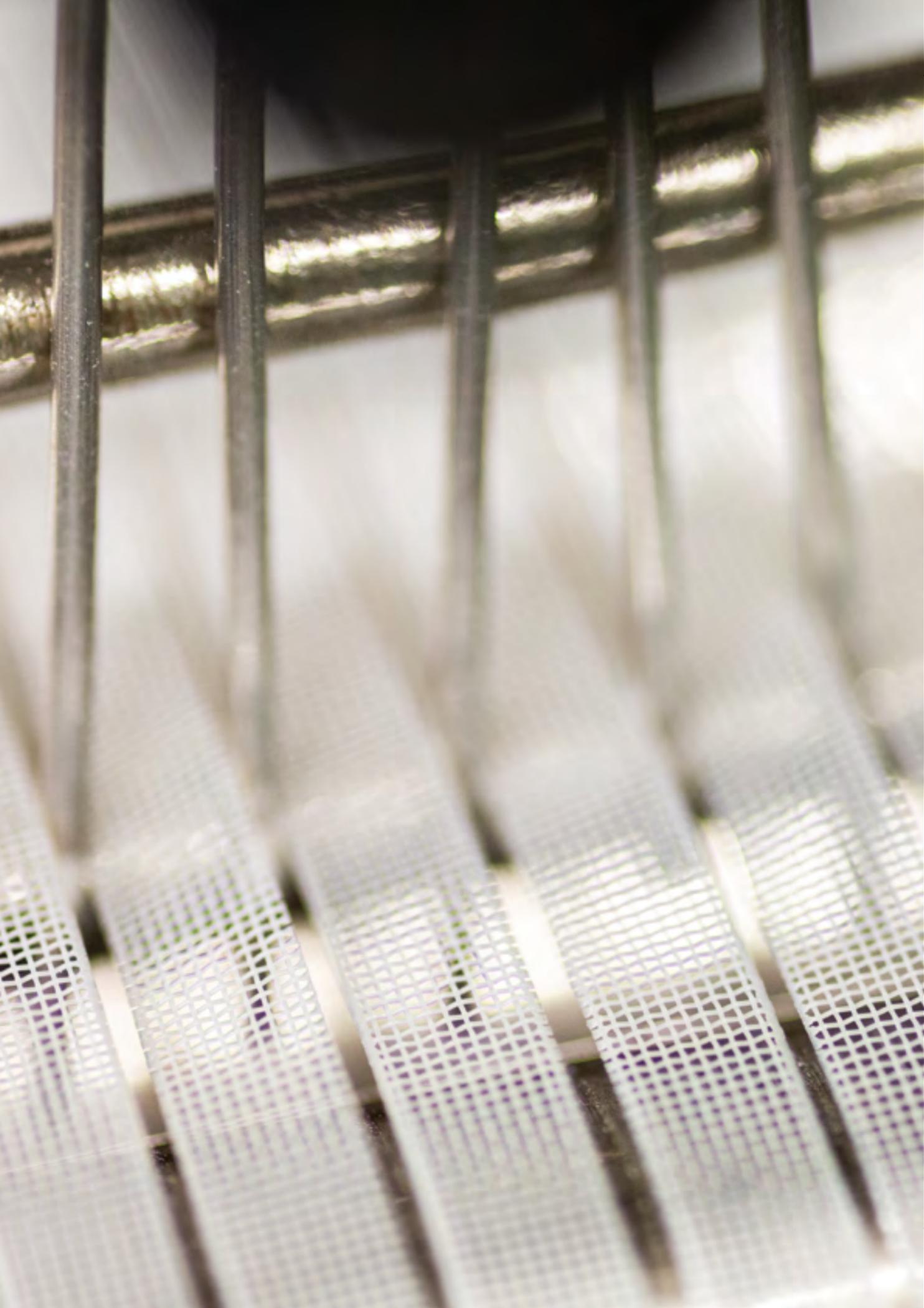

Ribbons

SAATCare fabrics can be heat-slit economically, or ultrasonic slit for the closest tolerances and the highest quality as needed.

Production Technology involved: Heat, Ultrasonics.

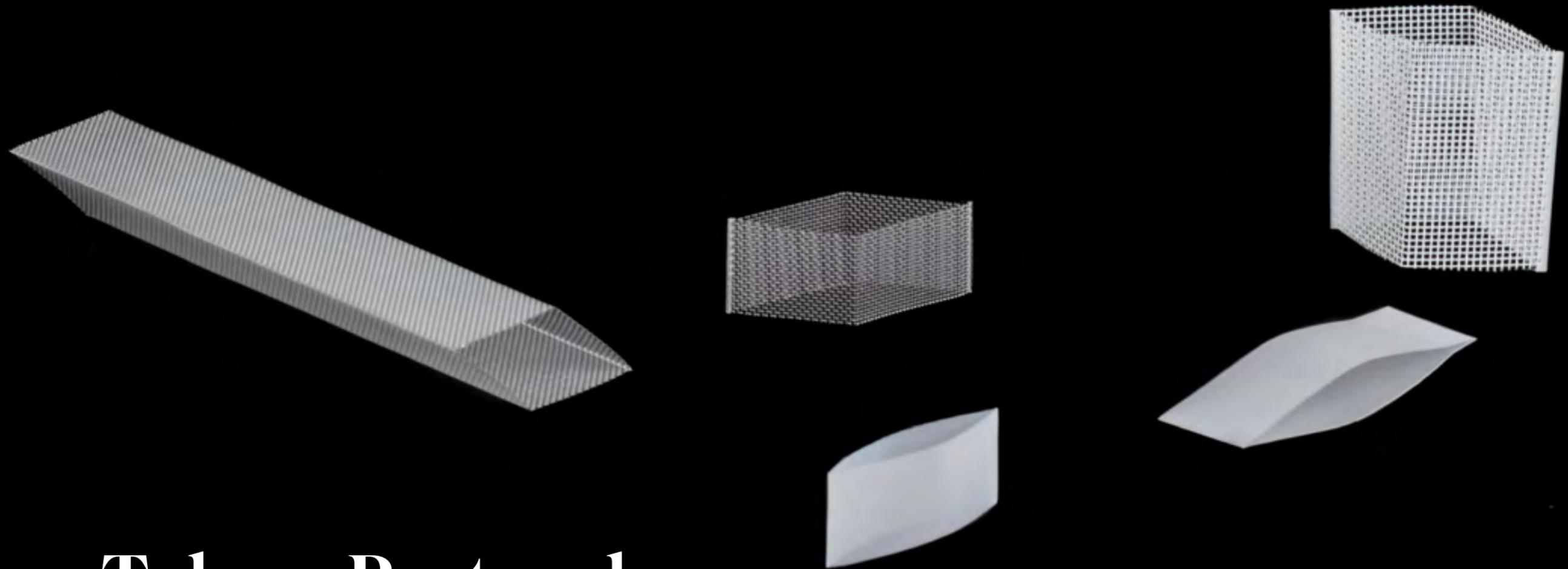

Tubes - Rectangles


For all applications requiring molded cylindrical filters, SAATI offers fabric tubes and rectangles with two open ends, die-cut or laser cut to length.

Many applications can accept the quality of a cold cut tube in view of its economic advantages. When a component must have one end sealed and one end open, SAATI is able to combine technologies in the same process and supply tubes with one end die-cut or laser cut and the other end ultrasonic sealed.

Ultrasonic technology is also applied for the fabrication of rectangular filters.

Production Technology Involved:

Die - Cut, Heat and Ultrasonic.


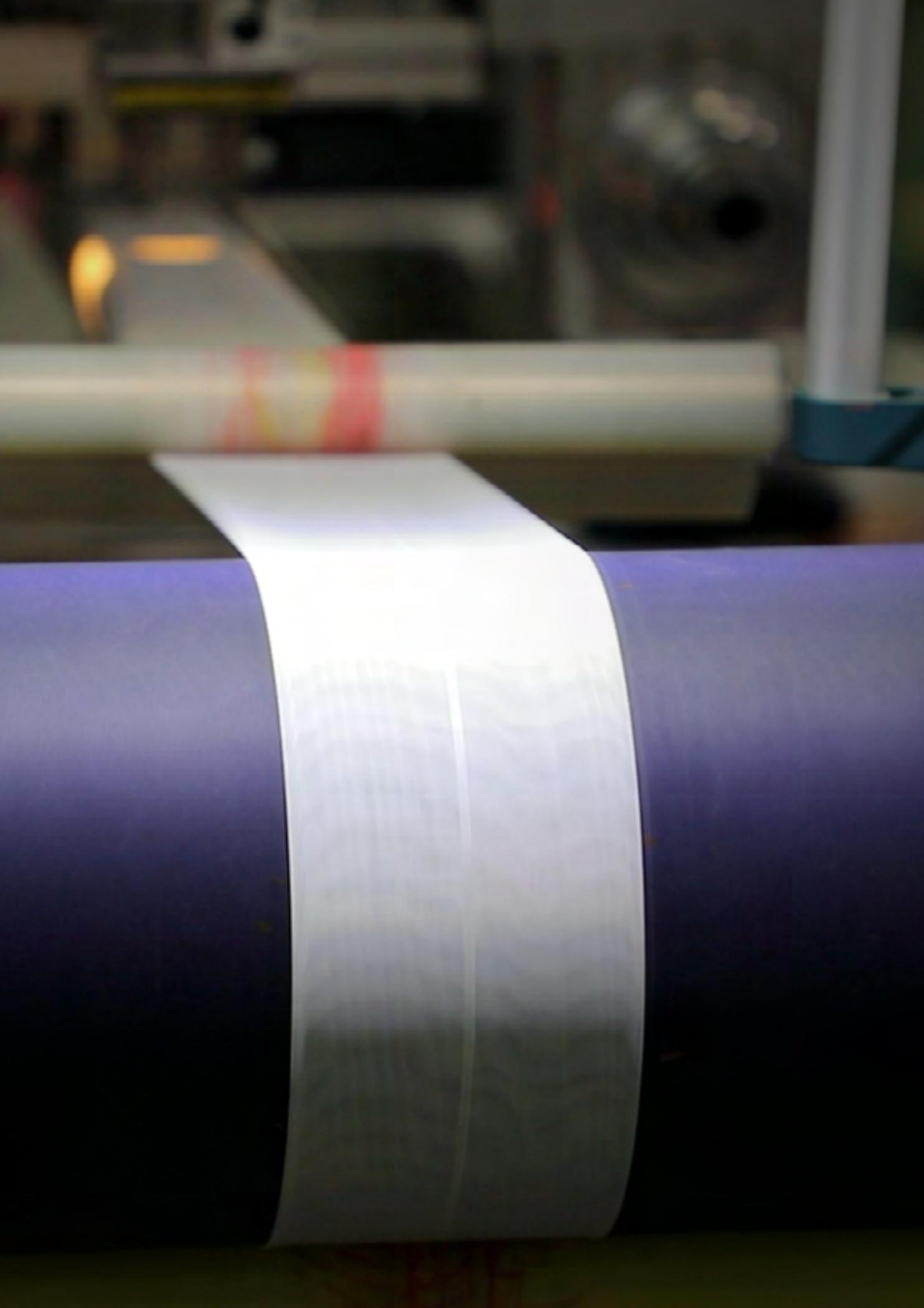


Continuous Tube/ Multi-Layer Products

Two layers of filter media (identical or different) are simultaneously slit with heated blades to form a continuous double-seamed tube.

Ultrasonic-welded tubular ribbons, although similar in construction to heat slit items, can be produced in a wide range of sizes, including very small ones. Two or more narrow layers can be attached using ultrasonic slitting. A fine filtration media can be supported or protected with a coarser one.


Production Technology involved: Heat, Ultrasonics.

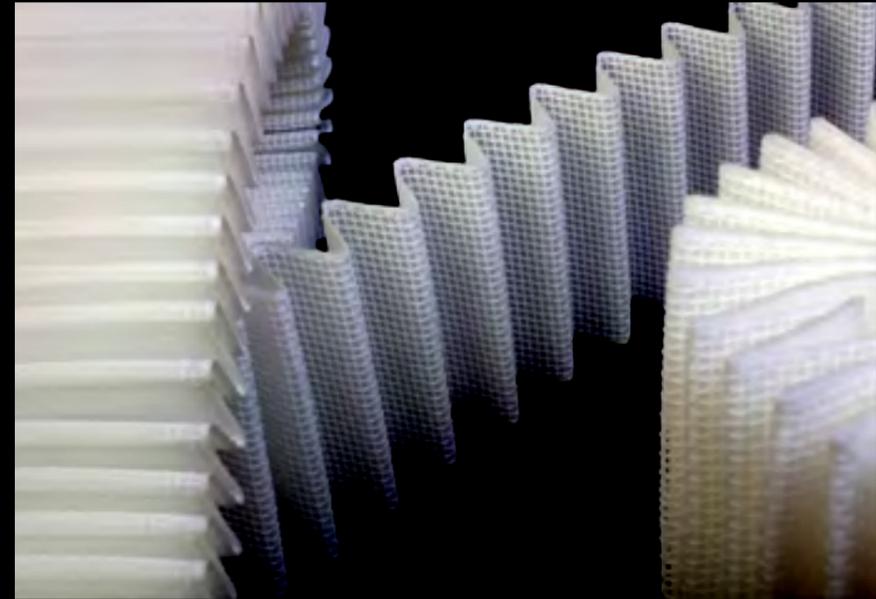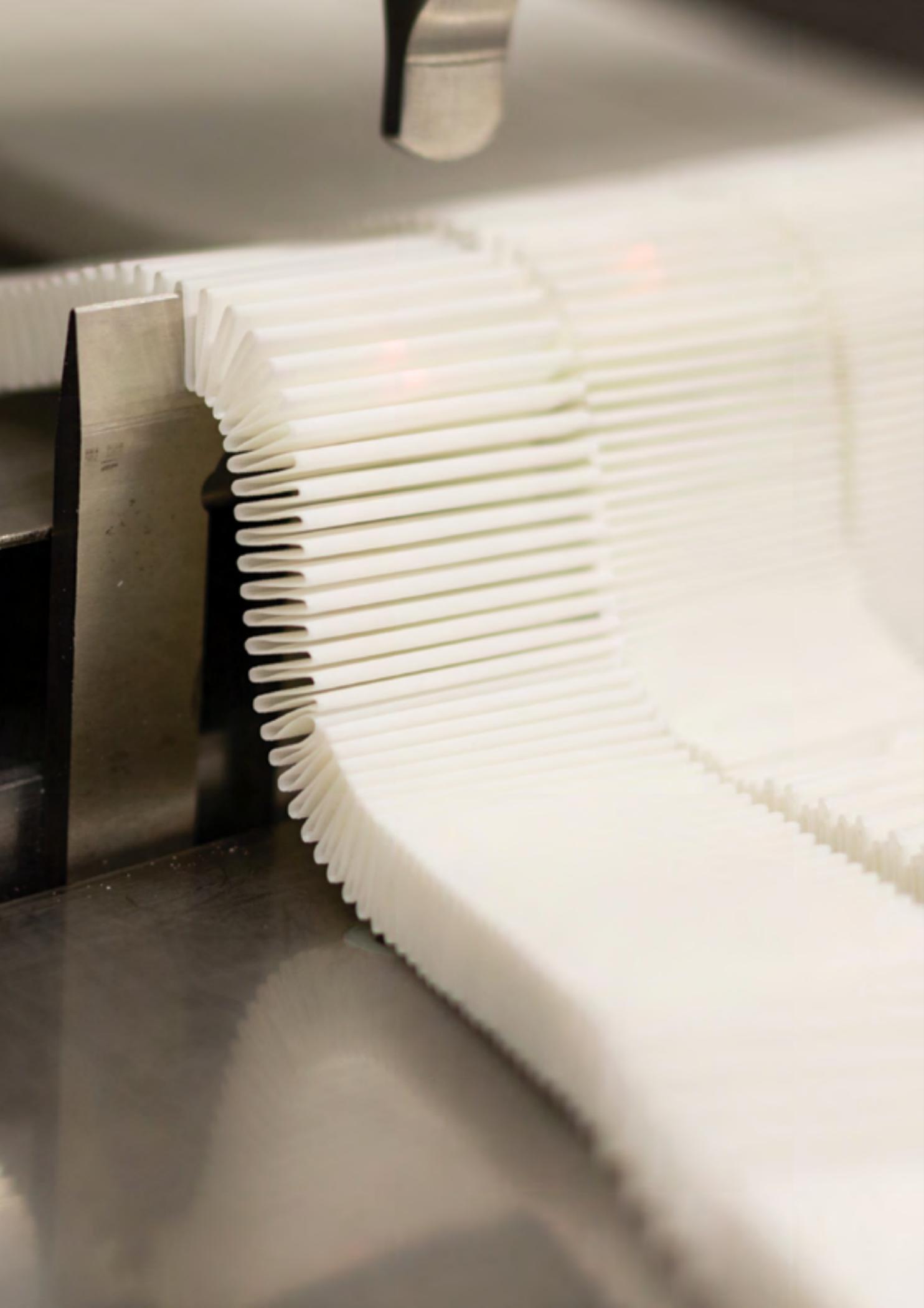

Pleated Components

Mono or multi-layer pleated components such as pack, ribbon, and cartridges can be manufactured for all applications requiring high filtration capacity in a narrow space.

Production technology involved: Heat, Ultrasonics.

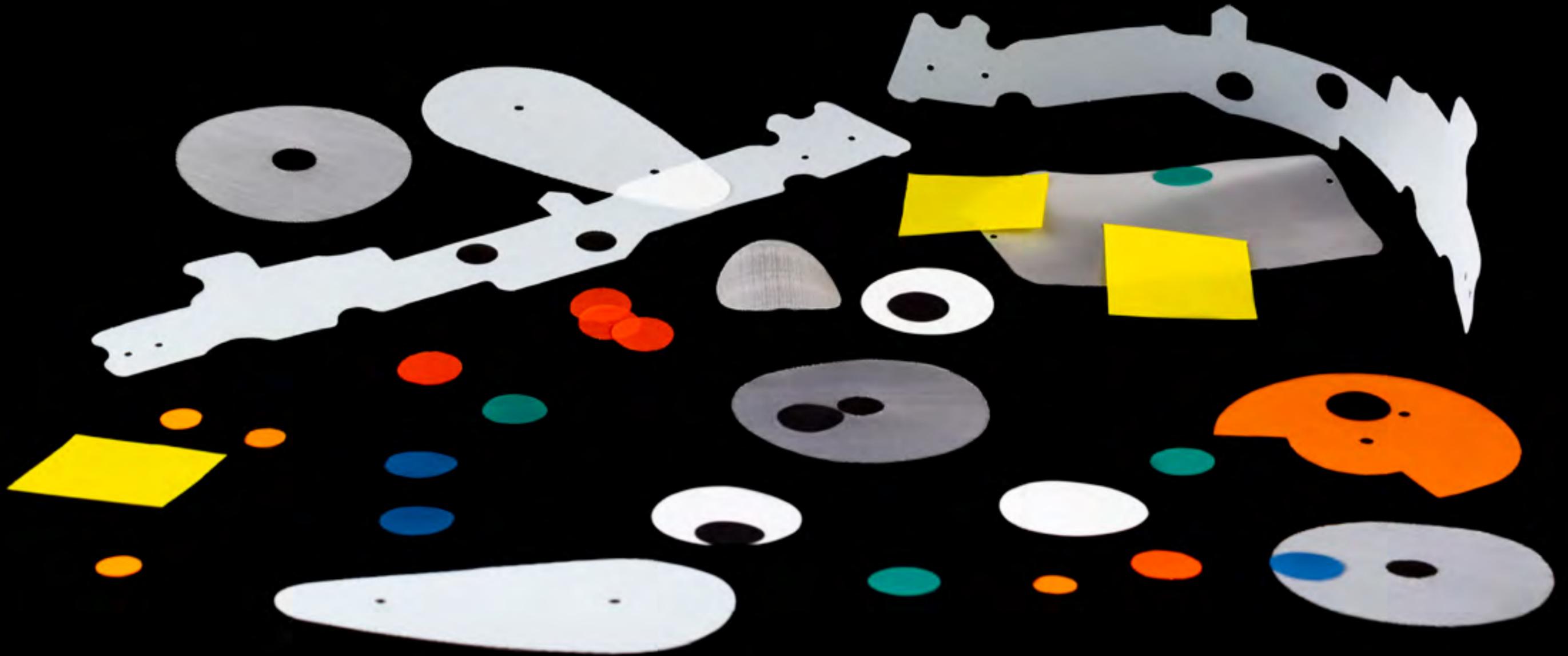

Custom Shapes

Die-cut or laser cut parts can be produced in almost any shape or size. One or more layers of fabric can be ultrasonic cut or sealed into virtually any shape using a PC controlled plotter, assuring a faithful reproduction of design.

Production Technology involved: Die-Cut, Ultrasonics, Laser.

Double Packaging

For all fabricated parts destined for use in medical applications, SAATI provides double packaging in order to ensure protection.

SAATicare Medical

High-Performance Features For Medical Applications

Precise pore sizes provide exceptionally selective filtration. Particulates and air bubbles of a specific size can be removed without affecting red blood cells. Materials with high open area provide high flow rates with minimal pressure loss. The uniform pressure drop provides consistent air pressure monitoring in spirometers. Smooth monofilament fibers and straight-through flow paths reduce the risk of hemolysis. Fabrics for medical applications can be specially finished to reduce pyrogens and contaminants.

-- SAATicare fabrics conform to USP Class VI

-- Fabrics are safe for gamma sterilization and autoclaving

SAATicare products designed for medical applications are exclusively medical grade.


Hemofiltration

Hemoperfusion, Hemodialysis, and Hemofiltration are methods of filtering the blood extra-corporeally to remove toxins or for blood purification in general. SAATicare woven meshes are used as protector filters in order to impede purifying elements (i.e. activated charcoal, resins, bicarbonate etc.) from flowing out of blood purifying cartridges.

Cardiosystem


During open heart surgery, heart and lung functions are temporarily replaced by medical devices (heart-lung machine) to guarantee the survival of the patient. Filters play an extremely important role in the heart-lung machine.



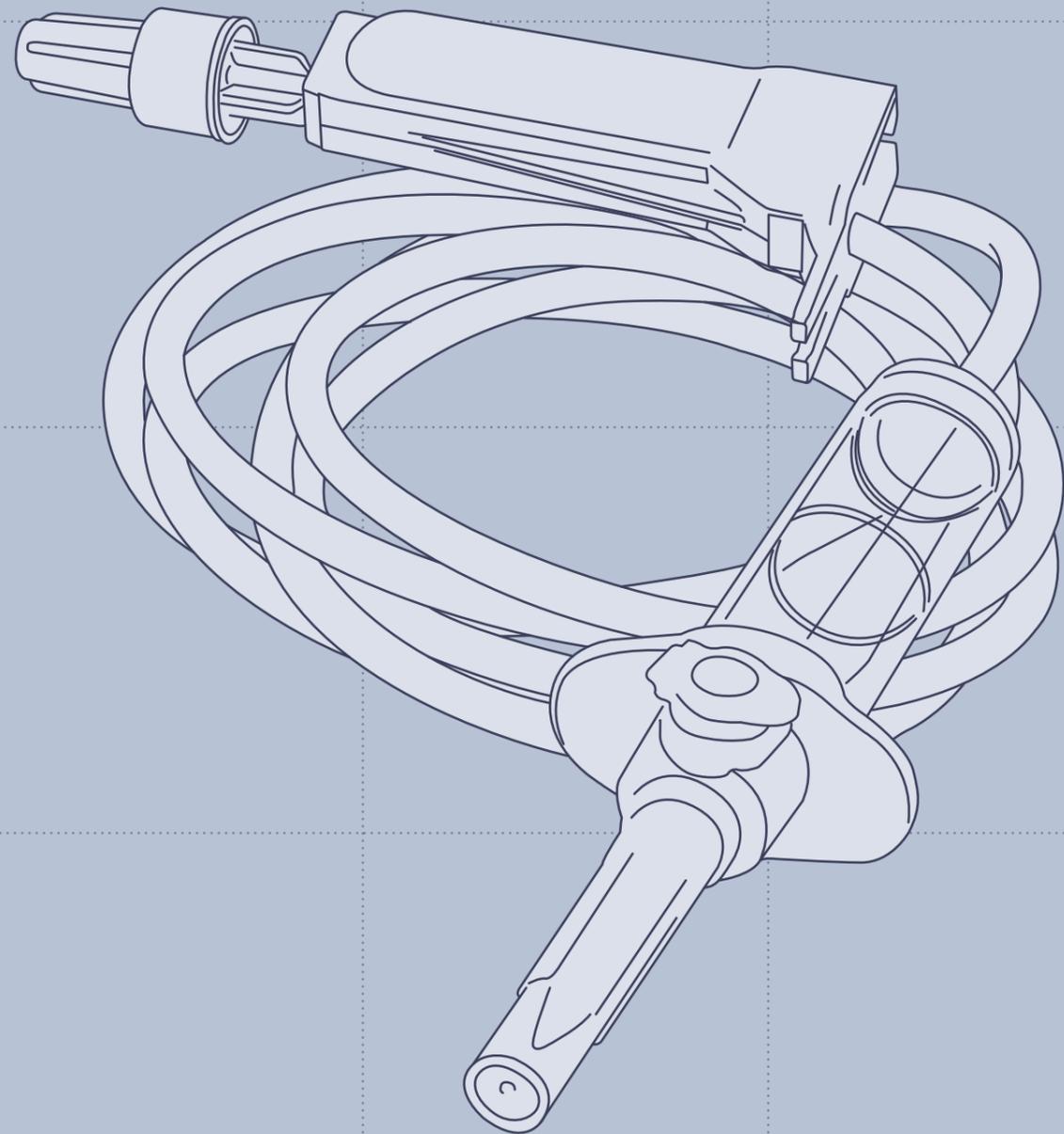


Transfusion Filters (Transfusion IV Sets)

SAATlcare products capture clots, micro and macro aggregates which could be present in the blood of the patient and reduce the risk of an embolism.

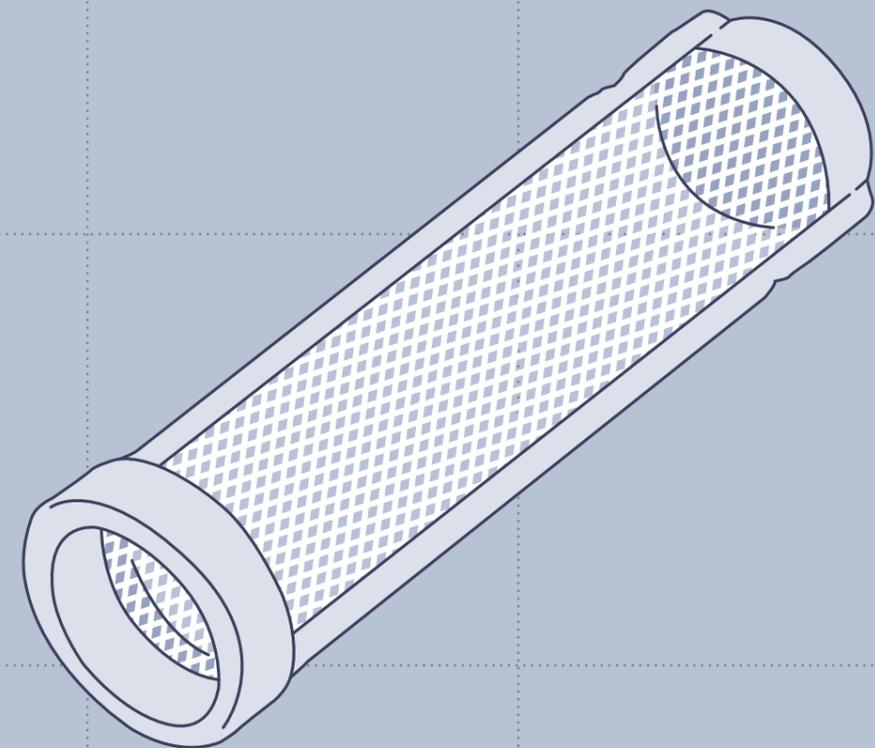

Infusion Filters

Possible contaminants originating from infusion bags may complicate the recovery of patients. SAATlcare monofilament woven meshes are designed to prevent possible blood poisoning caused by contaminants such as extraneous salts, microorganisms (protozoa, bacteria), dissolved impurities and other micro-particles.

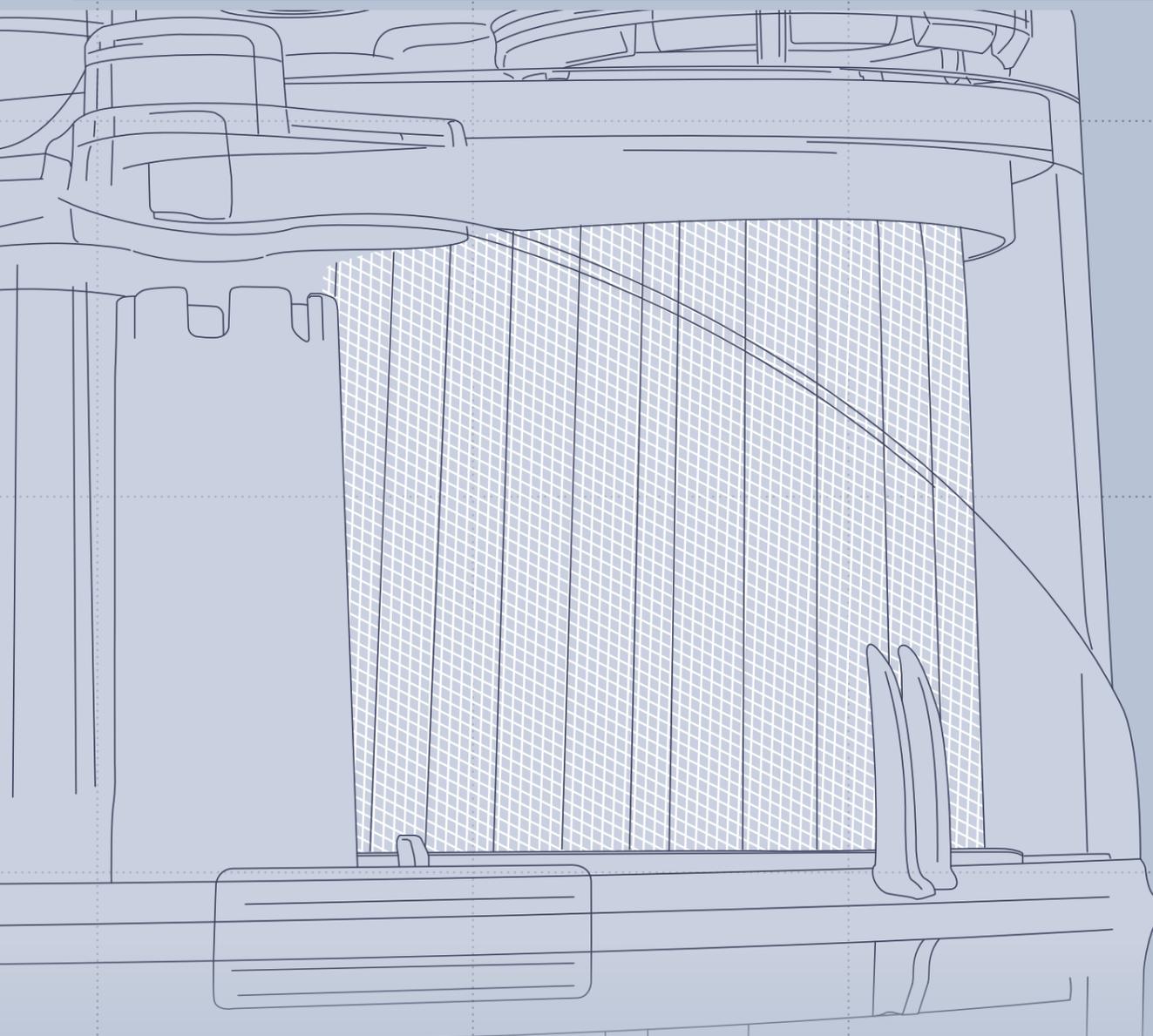



Cardiotomy Reservoir Filters

Possible foreign and undesirable particles such as thrombus, bone debris or air bubbles are harmful for the patient during surgery. This is the reason why a fine filtration is necessary. Furthermore, thanks to smooth monofilament yarn and straight-through flow paths SAATlcare products reduce the risk of blood cell damage.

Microaggregate Filters


The storage of whole blood is associated with formation of a considerable quantity of microaggregates between 25 μm and 40 μm . These microaggregates are mainly composed by degenerated leukocyte and thrombocyte fragments of red blood cells and other cells, lipids, lipoproteins and precipitated denaturated albumin. These microaggregates can damage pulmonary micro-circulation and for this reason SAATlcare meshes are used as filter media.

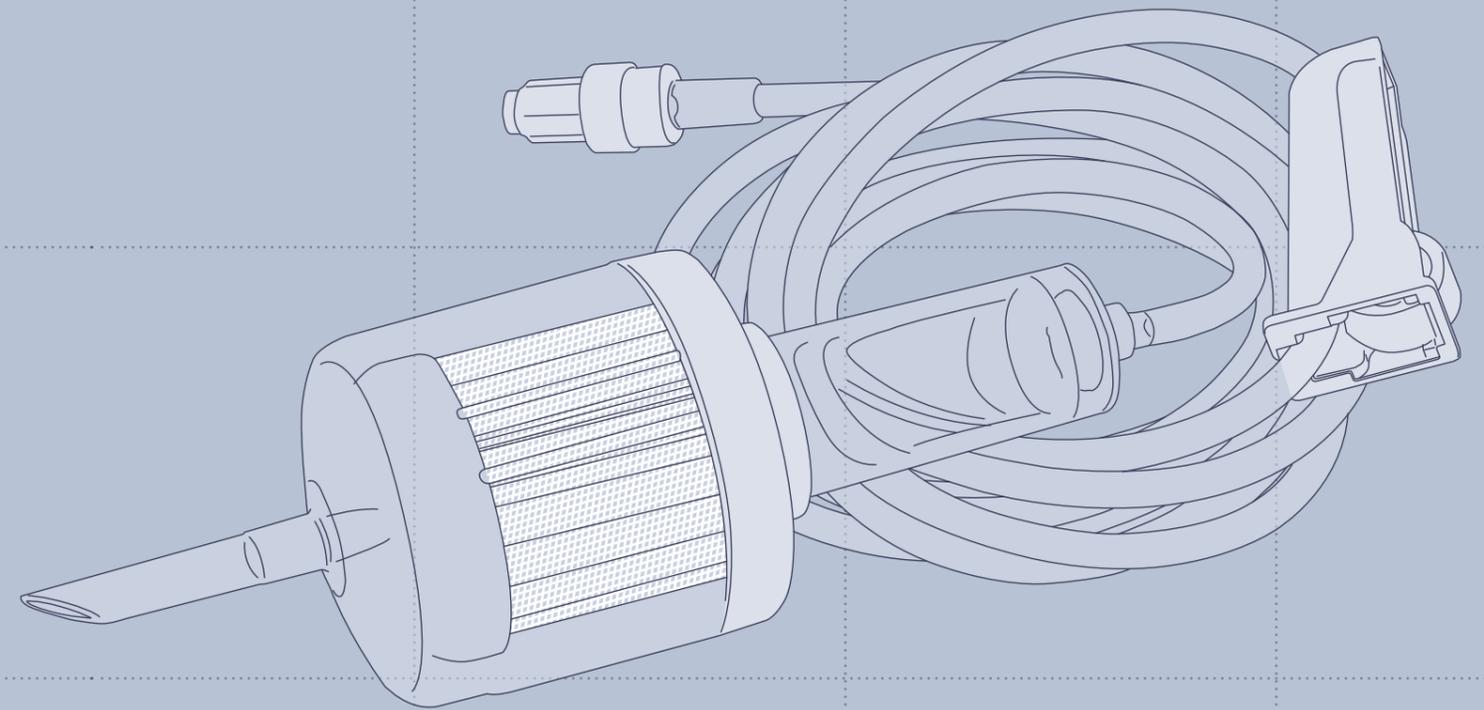

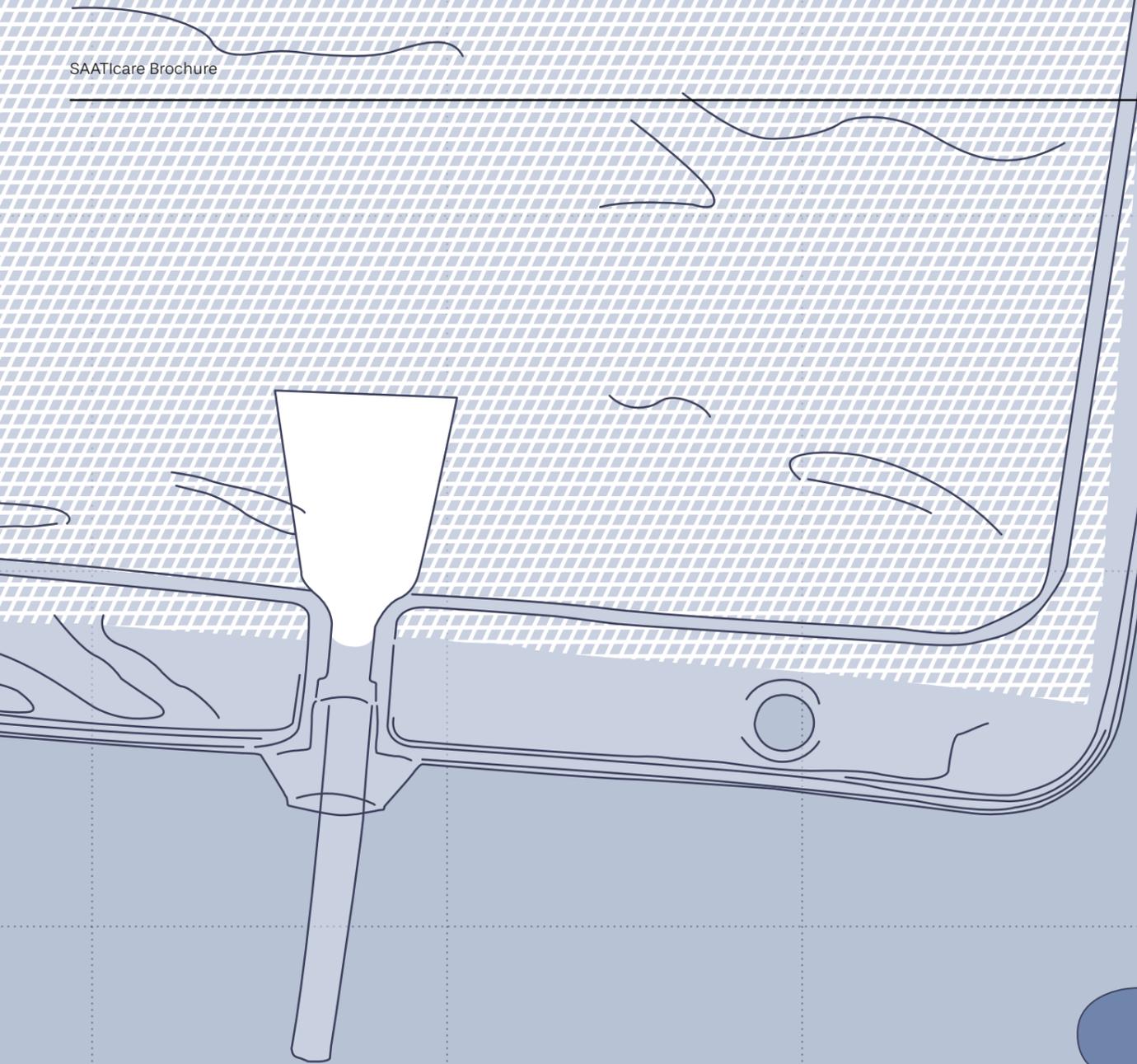

Blood Bags

The flow of blood must be filtered during transfusion to remove clots and small clumps of platelets and white blood cells that form during collection and storage.

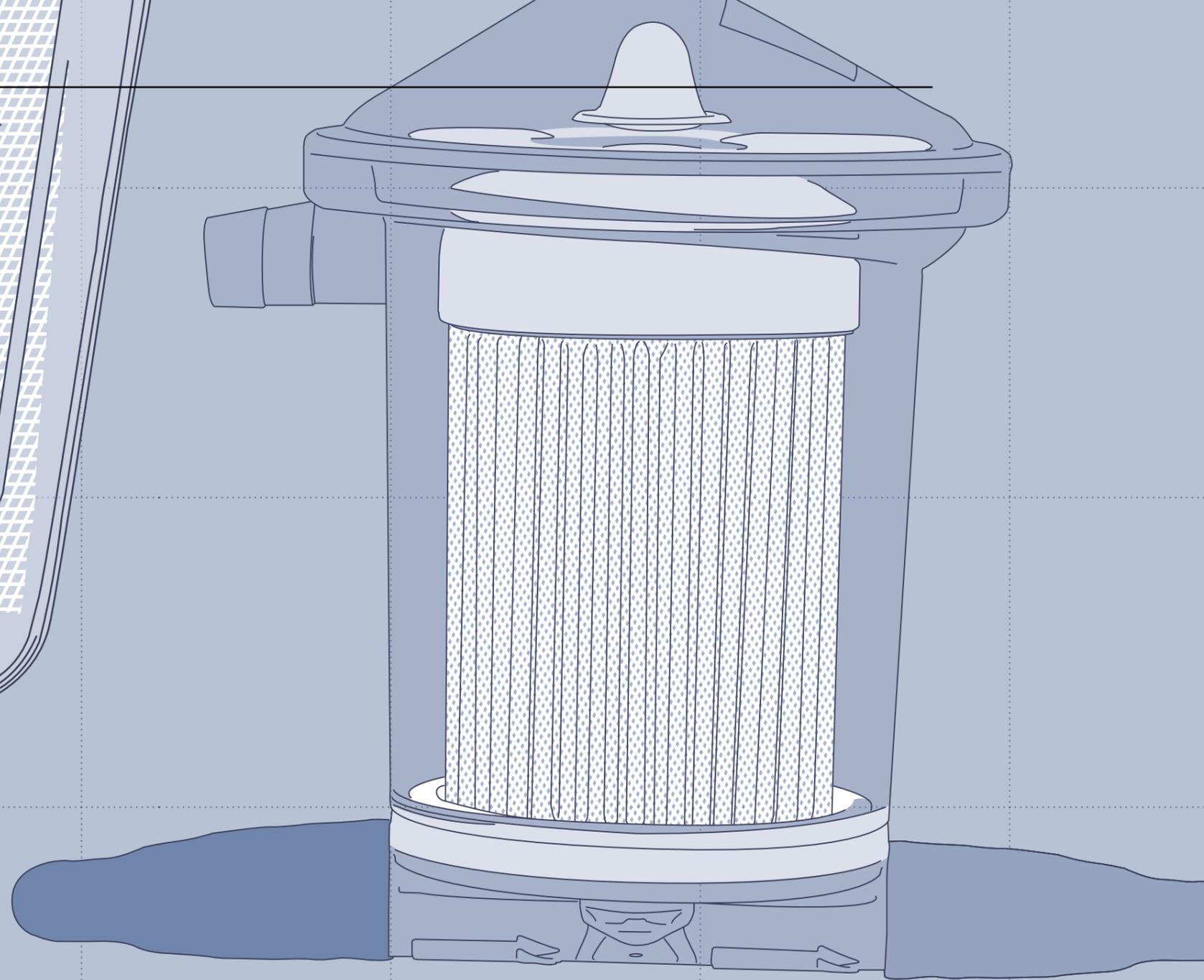


Arterial Filters

The aim of SAATicare mesh is to remove blockages in the perfusion circuit such as gas emboli, fat emboli, aggregates of platelets or red blood cells, air bubbles and other debris during the last filtration stage before the re-injection of the blood to the patient.

SAATicare Hyphyl™

A Surfactant Coating With Superb Hydrophilic Properties


SAATicare Hyphyl is a unique hydrophilic surfactant coating (an FDA-approved organic salt) that dramatically increases the amount of liquid transferred to lower layers of a test strip, as well as increasing the spreading area of a drop of liquid.

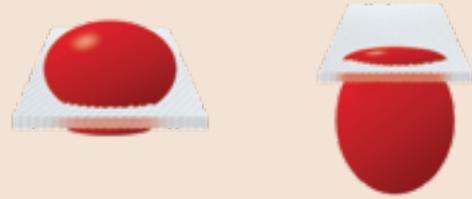

SAATI Hyphyl-treated Mesh makes it possible to perform testing with smaller samples, including blood samples




Low wicking observed in standard meshes

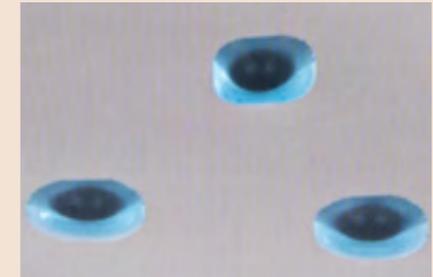


Excellent wicking observed in SAATicare Hyphyl-treated Meshes

Dynamic Contact Angle: Hyphyl-Treated vs Standard Fabric

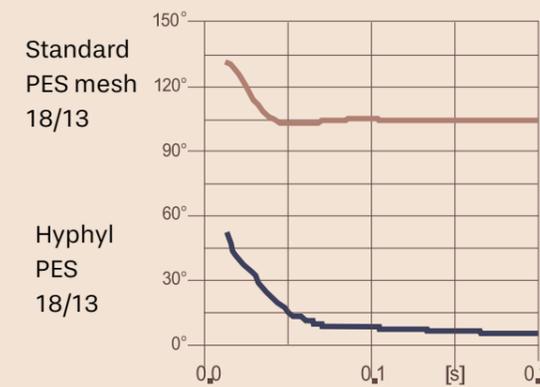

The lower contact angle of SAATicare Hyphyl indicates that it is more hydrophilic than standard mesh

Wet Area: Hyphyl-Treated vs Standard Fabric

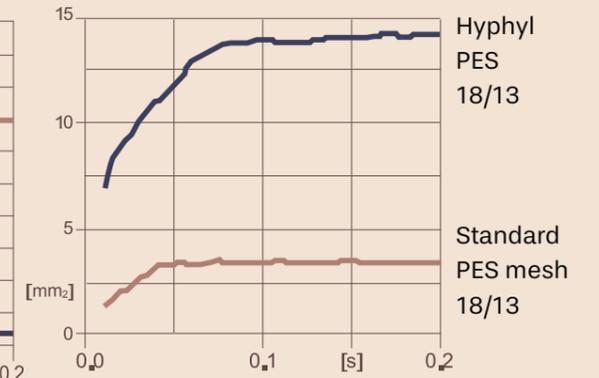

A 4-microliter drop of de-ionized water can wet an area almost five times greater using Hyphyl-treated fabrics


High Flow Fabrics



SAATlcare High Flow Fabrics for Arterial Line Blood Filters


SAATI has developed High Flow fabrics that achieve a superb higher open area compared to standard materials while maintaining the same mesh openings. For example, by using our SAATlcare High Flow fabric, the performance of arterial filter and other blood filters is increased up to 60% compared to standard fabrics. By guaranteeing a mesh opening even narrower 40 µm, such a high open area rate can dramatically improve the behaviour of the final filter, not just by minimizing pressure drop and priming volume, but even allowing a possible reduction in the filter's dimensions, with a considerable reduction of final cost.

Standard Fabric

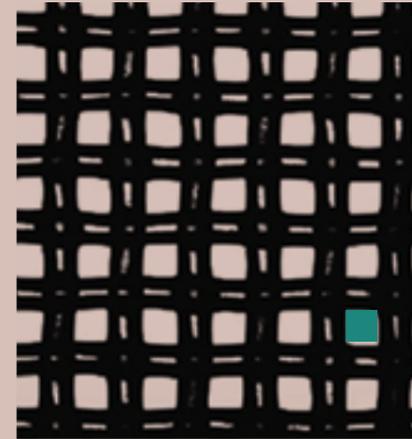


10x Magnification


SAATlcare High Flow Fabric

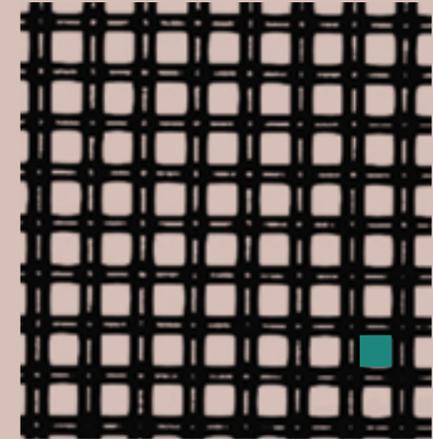

10x Magnification

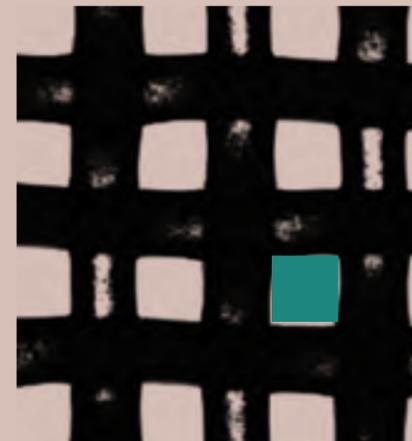

20x Magnification

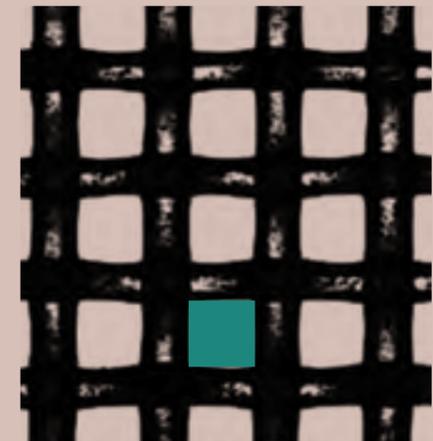

20x Magnification





SAATicare Diagnostic


High-Performance Screen Membranes For Critical Diagnostic Tests

SAATicare products are a complete solution for diagnostic test strips. They can be used as:

- Pre-Filter: Protective covering of the strip
- Support: Improving the reinforcement of the membrane and increasing its mechanical resistance thanks to the high tensile strength of the fabric. Use of fabric instead of non-woven material avoids the migration of the latter into the former as accurately as possible.
- Sample Pad: Thanks to SAATI fabrics treated with Hyphyl®, a special hydrophilic surfactant coating, the amount of liquid transferred to lower layers of test strip is dramatically increased. The precise pore size and weave provide high wicking rates for faster assay performance.

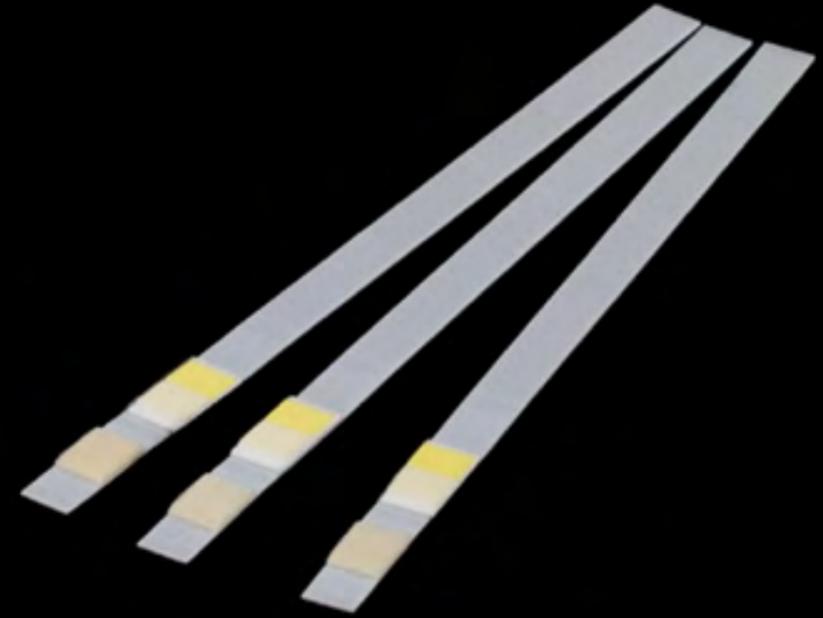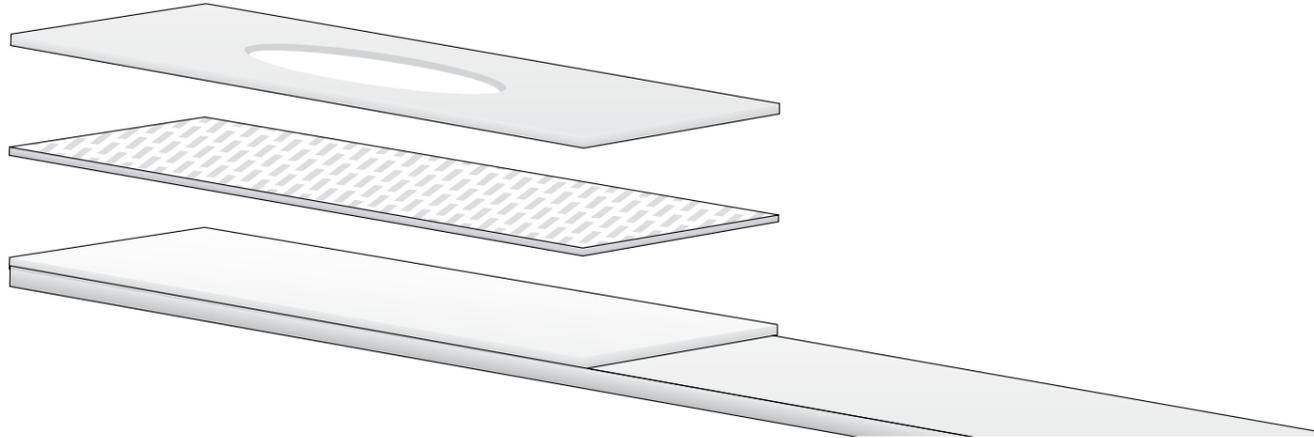

SAATicare fabrics enable diagnostic tests to be executed at their maximum performance level.

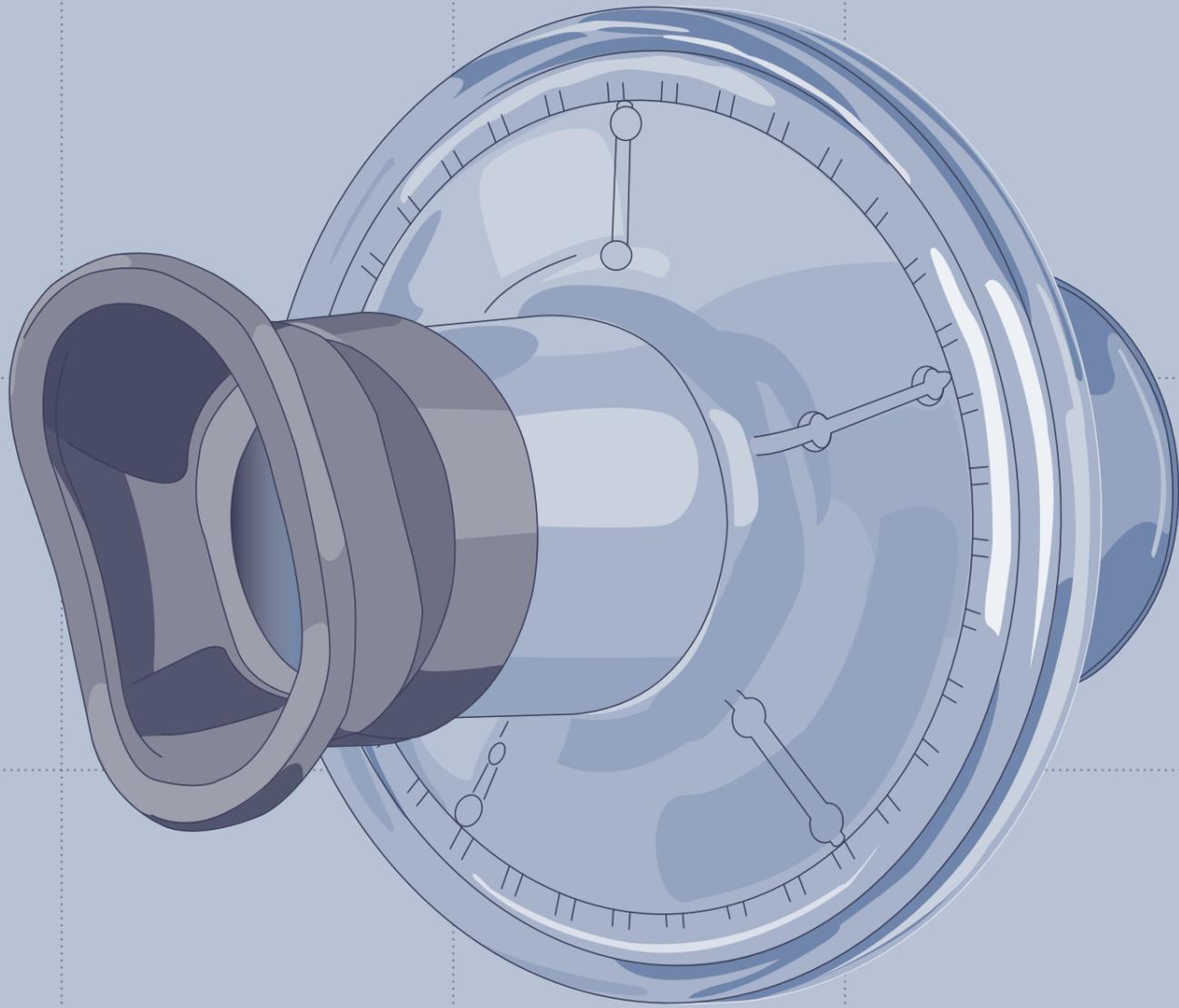


Spirometers


A spirometer is an apparatus for measuring the volume of air inhaled and exhaled by the lungs. The mesh is used as pressure transducer to measure pressure drop across the filter thanks to the consistency of air permeability.

Biopsy Bags

Biopsy bags are made of solvent-resistant material and are designed to process small biological specimens. Moreover they also maximize liquid exchange.

SAATlcare® biopsy bags can be produced with or without flaps.

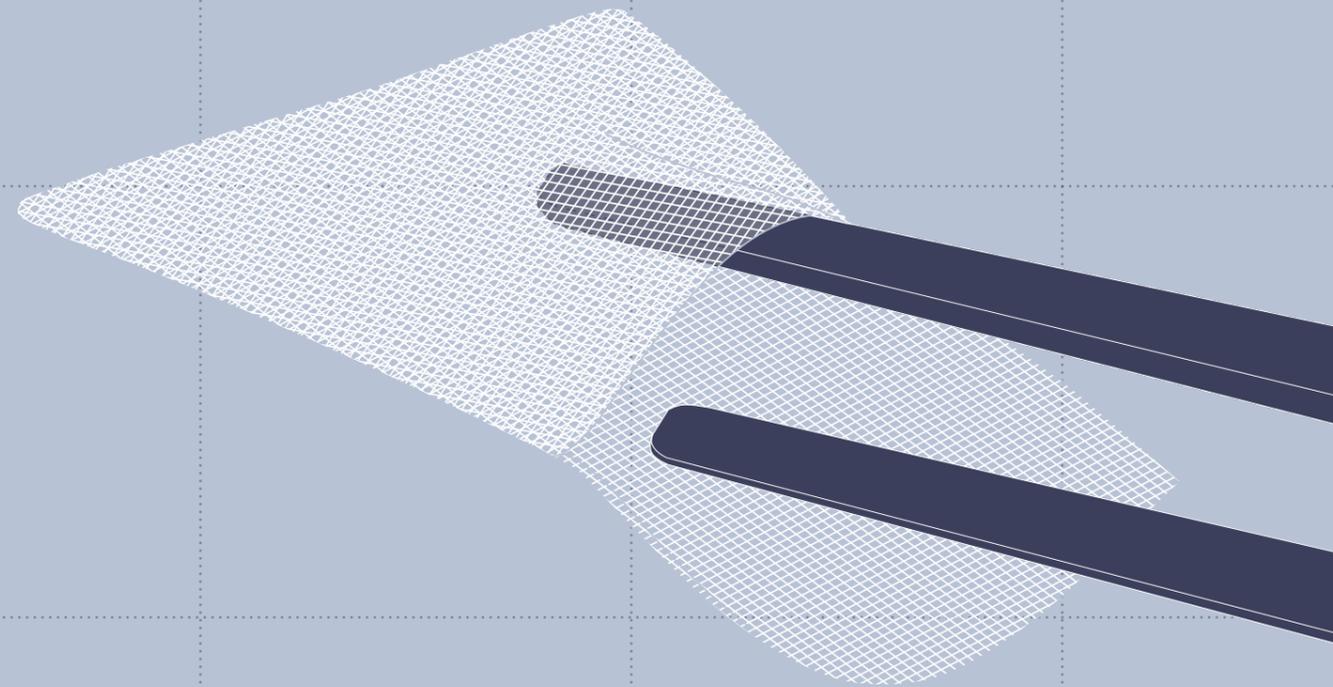


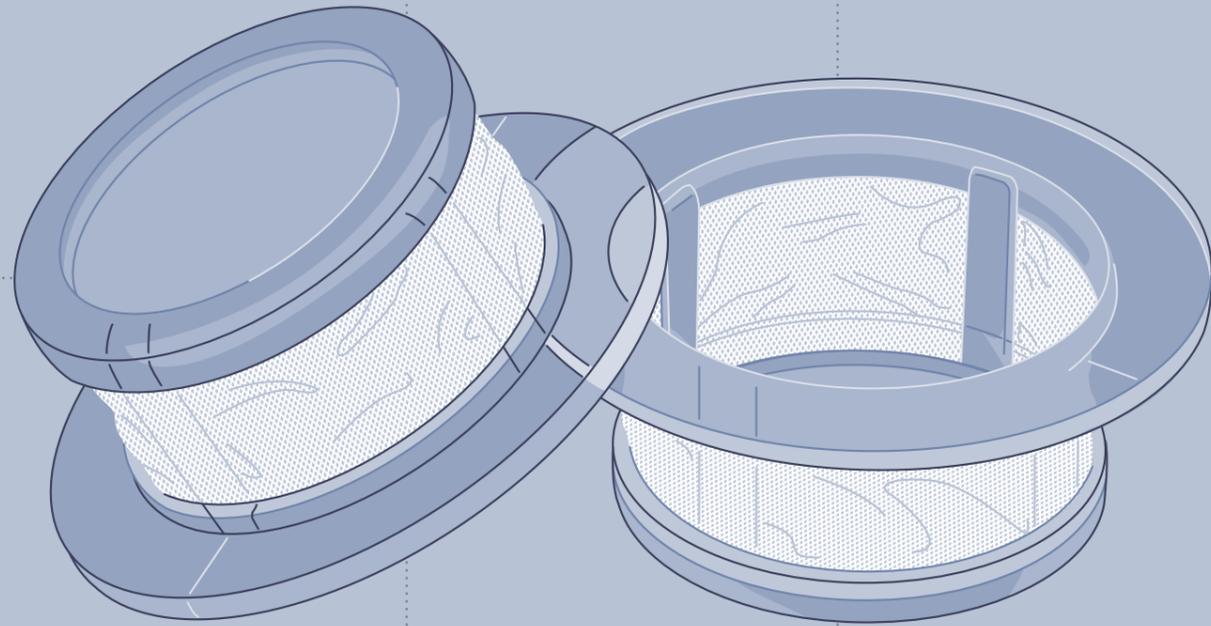


Molecular Sieves

Molded cups which use mesh with a specific pore size. It rapidly isolates primary cells in order to obtain a uniform single-cell suspension from tissues. SAATlcare mesh makes the process faster and easier compared to gauze filtration. SAATlcare's product range is characterized by a wide range of pore sizes. This makes a perfect media for sieves used during preparation of flow cytometry samples, specimens for primary cell cultures/immunogens and freezing stocks. SAATlcare products filter also agglutinative proteins produced in inactivation serum.

Other Applications

SAATlcare products are also used for fluid handling, sample preparation, reagent delivery, and as a **carrier for membranes**.

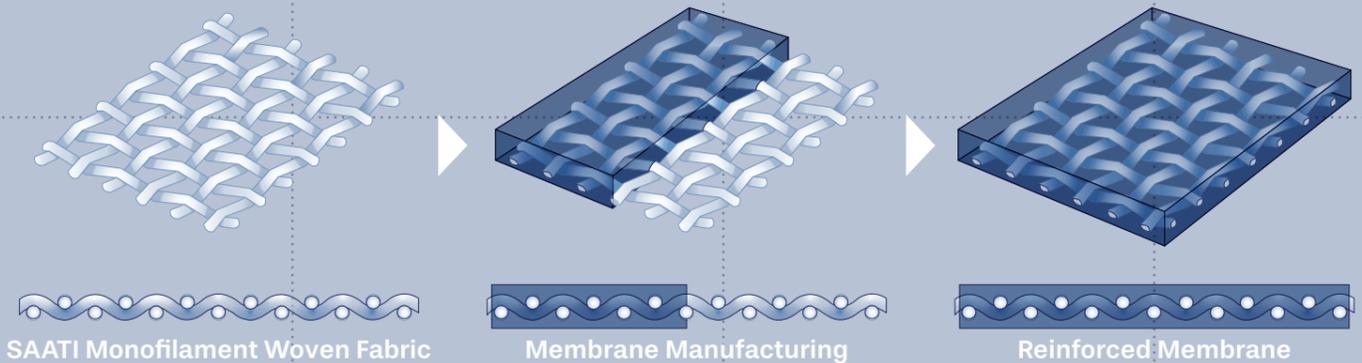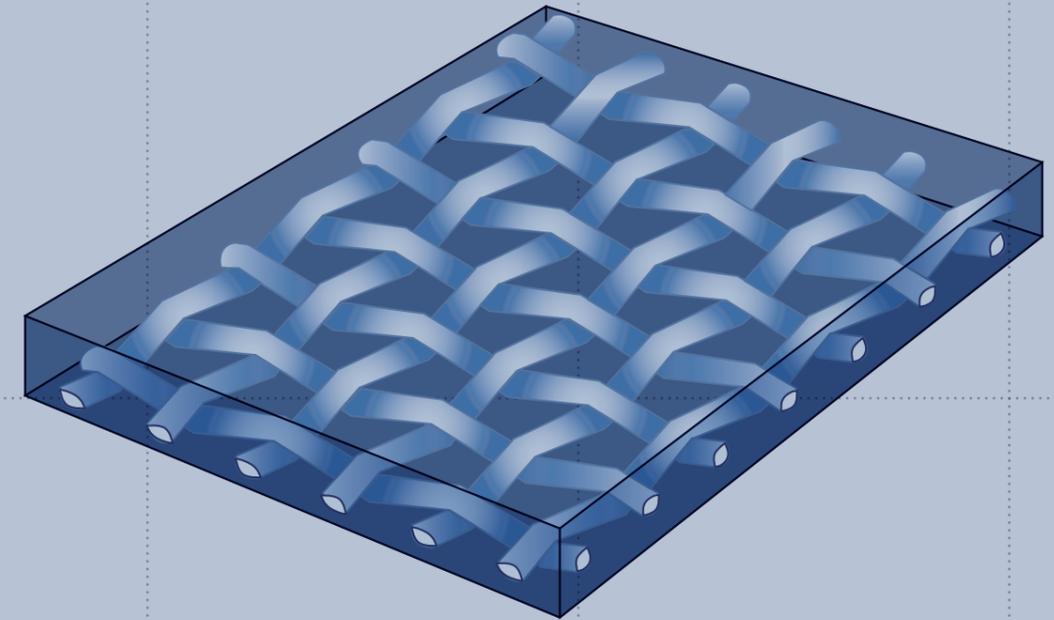

Testing Compliance

SAATI medical grade fabrics can be tested and certified according to the following regulations:





USP Class VI Test USP 24/NF 19, 2000

Biocompatibility of plastics. Including Systemic Toxicity, Intracutaneous Toxicity and Implantation Tests.

Hemolysis Test

DIN 58970

Detecting compounds which may damage red blood cells.

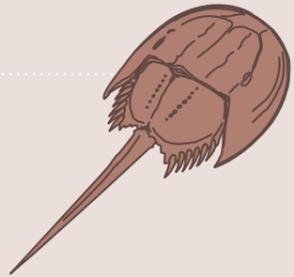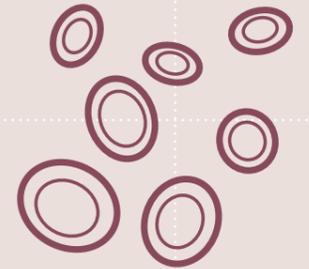

LAL Test

USP 24/NF. 19, 2000

Limit: 0.125 EU/ml



Pyrogen Endotoxin Test.

Detects endotoxins released after the death of bacteria.

Cytotoxicity Test

ISO 10993 Part 5

Investigates whether or not the leachables from the material may cause cytotoxicity (death of cell).

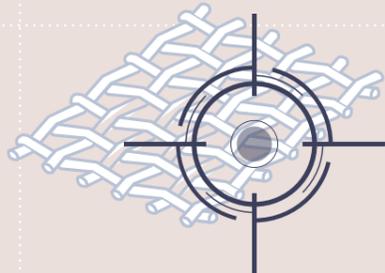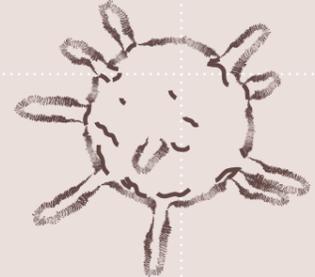


Extractables

21CFR177.1630 21CFR177.1500

Amount of extractables from a filter during its intended use.

Additional Tests

Additional tests are available on request.

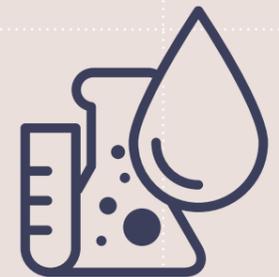

A doctor in a white coat with a stethoscope around his neck is smiling and looking at a young boy lying in a hospital bed. The boy is also smiling and looking up at the doctor. The scene is set in a hospital room with blue walls and a window with blue curtains. Medical equipment is visible in the background.

Sustainability Approach

We Work To Improve Everyday Life

We understand the importance of operating a sustainable business. First of all, in terms of the products we manufacture. At SAATI, we develop and make available technological products that feature in people's everyday routines, all over the world, helping to make their lives healthier, safer, or simply easier.

A family consisting of a man, a woman, and two children are walking away from the camera in a grassy field. The scene is bathed in the warm, golden light of a sunset. The woman is carrying a wicker basket. The background shows a line of trees and a bright sky.

Life Is Precious. We Work To Protect It

All of our products, from fabrics, through components to chemicals, recognize that life is precious. At SAATI, we manufacture technologies for life: filter fabrics for blood transfusion sets; aramid fabrics in bullet proof jackets; functional fabrics and chemicals to manufacture solar cells, mobile phones and tablets that improve our lives.

SAATI S.p.A. -**World Headquarters**

Via Milano 14
22070 Appiano Gentile Italy
Phone: +39 031 9711

SAATI France S.A.S

74 Route de Bapaume - Sailly-Saillisel
80360 - France
Phone: +33 3 22 85 77 00

SAATI Serigrafia Iberica S.A.U.

Pl. "El Mijares" c/Industria 13 - 12550
Almazora (Castellon) - Spain
Phone: +34 964550688

SAATI Americas Corp.

201 Fairview Street Ext.
Fountain Inn, South Carolina 29644 USA
Phone: +1 (800) 431-2200

SAATI Technical Fabric (Tianjin) CO. Ltd

Cross of Saida 2nd Branch Road, Saida Century
Avenue, Xiqing Ec.Dev. Area - Tianjin - China
Phone: +86 22 23960843

OOO SAATI Russia

Shvetsova Street, 23 house, 198095,
St. Petersburg, Russia
Phone: +79 062788343

SAATI Korea Ltd

SK Ventium 101-601 - Dangjeong-dong, Gosan-ro
166, Gunpo-si, Gyeonggi-do, South Korea, Zip 15850
Phone: + 82 31 429 9337

SAATI**Deutschland GmbH**

Ostring 22 - D-46348 Raesfeld - Germany
Phone: +49 2865 95800

SPT - Sales + Marketing GmbH

Kurpfalzring 100a - 69123 Heidelberg - Germany
Phone: +49 62 21 | 77 876-27

CST Colour Scanner Technology GmbH

Koenigsberger Str. 117 - 47809 Krefeld - Germany
Phone: +49-2151-1592260

SAATI Advanced Chemicals LLC

201 Fairview Street Ext.
Fountain Inn, South Carolina 29644 USA
Phone: +1 (218) 628-2217

WuXi TianYi Precision Fabrics Co., Ltd.

No. 28, Xigang West Road, Donghutang,
Donggang Town, Xishan District, Wuxi City - China
Phone: +86 510-88791064